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
Model: RH-700122

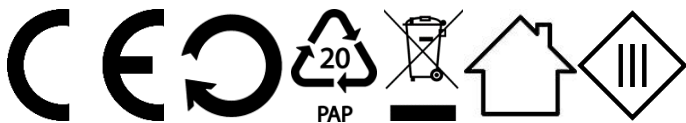SPIS TREŚCI

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Model: RH-700122

Table of Contents

1. INTRODUCTION	12
2. INTENDED USE	12
3. TECHNICAL DATA	12
4. PRINCIPLES OF SAFE USE	12
5. EXPLANATIONS OF THE SYMBOLS	14
6. CONTENTS OF THE SET	15
7. CONSTRUCTION	16
8. USE	16
9. CLEANING AND MAINTENANCE	18
10. TROUBLESHOOTING	18
11. STORAGE AND TRANSPORT	19
12. WARRANTY AND SERVICE	19
13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS	19
14. MANUFACTURER	19

1. INTRODUCTION

Thank you for choosing our product. We are confident that its use will be a great pleasure for you. Please read the user manual carefully before using the product, in particular the instructions on safe use. Please keep this manual in order to be able to refer to it while using the product. The manual is also available in electronic form on the manufacturer's website: <https://mpotech.eu/qr/hoffen>

2. INTENDED USE

The device is designed for massaging lower back, back , neck and thighs. This device is intended for private use only in a house or car, and may not be used for commercial purposes. Any other use is improper and dangerous. The manufacturer is not liable for injuries resulting from the use of the device in a manner inconsistent with its intended use. The device is not a medical equipment.

3. TECHNICAL DATA

Name	Pillow shiatsu massager
Model	RH-700122
Lot number	202110
Power	20 W
Power supply	12 V = 2 A
Power adapter	Input: 100-240 V ~ 50/60 Hz, 0,8 A Output: 12 V = 2 A, 24 W
Additional	2 heads with 4 massaging points, a rubber band to stick to a headrest, car power cord

4. PRINCIPLES OF SAFE USE

Read these guidelines carefully. Not following them may be dangerous or illegal.

–IMPORTANT SAFETY INSTRUCTIONS–

- This appliance may be used by children from 8 years of age and older and persons with limited physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised or have been instructed in the safe use of the appliance and understand the hazards involved.

- Children must not play with this device.
- Children must not clean or maintain the appliance without supervision.
- Children under the age of three must not use this device due to its lack of response to overheating.
- The product can only be used with accessories supplied with the device.
- If the power cable or car power cable is damaged, it must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid a hazard.
- The surfaces of the unit might heat during work. People whose skin is not sensitive to heat should exercise caution when using the device.

—DO NOT RISK—

The device is designed exclusively for indoor use in a house, an apartment or a car. Not suitable for outdoor use (terraces, balconies, etc.). Do not use this device in locations where it is prohibited to use electrical equipment, or where it may cause harmful interference or other hazards. Always connect the device to a power supply that matches the data on the rating plate.

Do not switch on or use the device close to flammable materials. Do not place any flammable objects on the device - this may result in damage to the device or even fire. Do not lift or carry the device while it is switched on.

—SURROUNDINGS—

Do not use or store the device in dusty, damp or hot conditions. Keep children under 8 years of age and pets away from the unit when it is in operation, as the surfaces of the unit become very hot during operation - there is a risk of burns. The power cord may cause tripping, tipping over or other injuries. Do not leave working device unattended.

—WATER AND OTHER LIQUIDS—

Do not expose the device to water or other liquids. It can influence work and the life span of electronic components. Avoid working in an environment with high humidity - steamed kitchens, bathrooms or drops

EN © 2022 HOFFEN. All rights reserved.

Pillow shiatsu massager

of water may cause moisture to appear in the device, which may result in damage to the device and in extreme cases of electric shock. Do not use the device when your hands are wet or damp.

–QUALIFIED SERVICE–

Only qualified personnel or an authorized service centre may install or repair this product. Repairing by an unauthorized or unqualified service may cause device damage and warranty loss. The manufacturer is not responsible for problems caused by unauthorized modifications to the equipment made by the user.

-CONTRAINDICATIONS-

Do not use the device without prior medical consultation in the following cases: pregnancy, blood circulation disorders, varicose veins, open wounds, bruises, skin splits, phlebitis, you have a pacemaker or joint prosthesis, you are using medical therapy.

5. EXPLANATIONS OF THE SYMBOLS	
	Important information on the use of the device.
	WARNING! Please read the warnings carefully. Failure to do so may result in damage to the device and, in extreme cases, in loss of health or life.
	Symbol for cardboard (packaging material).
	The mark placed on the product indicates that the product complies with the basic requirements and other provisions of the so-called "Directive of the 'New Approach' Directive of the European Union.
	The packaging of the device can be recycled in accordance with local regulations.
	The packaging elements should be disposed of in appropriate containers.
	Marking of the container into which the cardboard - paper should be disposed.

	Indoor use only.
	Symbol of Class III equipment. A Class III appliance is supplied from a separated/safety extra-low voltage (SELV) power source, in such system the voltage cannot exceed ELV under given conditions.
 	The device is marked with a crossed-out rubbish bin, in accordance with the European Directive 2012/19/WE on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.

6. CONTENTS OF THE SET

Included in the set:

- Pillow shiatsu massager
- AC Power Supply
- Car power cord
- User Manual

 *If any of the above items is missing or damaged, please contact your dealer.*

 *We advise to keep the container in case of any complaints. Keep the packaging out of reach of children.*

You print a version of the user manual with the large font size or download it from the web address below:

<https://mptech.eu/qr/hoffen>

7. CONSTRUCTION

- | | |
|---------------------------------|--------------------|
| 1. Case | 5. Massaging heads |
| 2. Rubber band for the headrest | 6. Car power cord |
| 3. On / Off | 7. AC Power Supply |
| 4. Connection hose | |

8. USE

 *Ensure that the cable is routed so that there is no possibility of damage, cutting, pushing or pulling that could cause damage to the equipment or injury to you or anyone else.*

 *If you notice that the device cover, the mains adapter or its cable or the car power cord is damaged, do not switch the device on.*

 *Be careful not to let water get inside the device! This is an electrical appliance and water may cause a short circuit that may damage the appliance or cause electric shock.*

BEFORE USE

Remove the unit from the packaging, unwind the AC adapter cable or car power cord and connect it to the connection cable on the back of the unit. Connect the power supply plug to the mains socket or the car plug to the cigarette lighter socket in the car.

 *Using the device without the cover in place is dangerous and may lead to injury. Use the device only with the case on.*


SWITCHING ON/OFF AND DEVICE OPERATION

To switch on the device press the on/off button, the device will start with the backlighting of the massage heads on. To switch on the device press the on/off button, the device will start with the backlighting of the massage heads on. Pressing the button 4th time will switch the device off.

In addition, the device switches off automatically after 15 minutes of operation.

When you have finished using the device, switch it off and disconnect the AC adapter or car power cord from the power source and the device itself.

MASSAGE

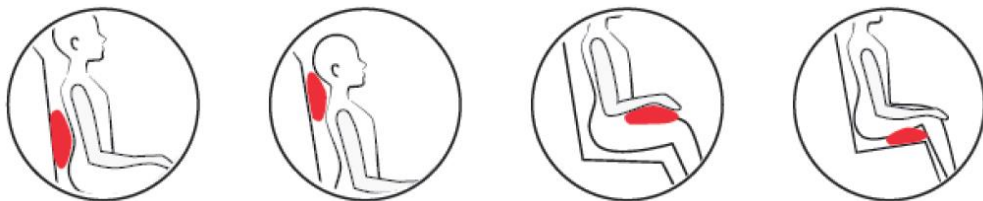Massage position

1. Depending on the body part being massaged, place the device on the back of a high chair, on your thighs or under your thighs. In the car, it can be attached to the headrest of the car seat using the elastic band on the case.
2. Sit in a straight position and make sure the device is positioned correctly.
3. Turn on the device to start the massage.

Tips on the massage:

- Do not massage one place at a time for too long a time, to avoid skin irritation and redness.
- In case of skin irritation stop massaging at once.
- Do not massage areas with irritated or damaged skin.
- Do not use the device around the area of eyes or other sensitive body parts.
- In case of any discomfort stop the massage at once.
- Do not use the device as a help in therapy or instead of medical treatment. The device is not a medical equipment.

9. CLEANING AND MAINTENANCE

 **WARNING!** The device must not be submerged in water or washed under running water!

Make sure your hands are dry before cleaning the unit. The device has to cool down and be disconnected from the power supply.

 *The device requires no maintenance.*

CLEANING

The case is fit for hand wash only. Unzip the zipper and of the case and take it off from the device to wash it. Wash in warm water using delicate detergents (i.e. watered-down dishwashing liquid). Do not bleach, do not clean with chemicals. Once washed do not tumble dry and do not iron.

10. TROUBLESHOOTING

In case of any problems with the device, please refer to the following tips.

The device does not switch on	Check that the power supply plug or its cable or the vehicle's power supply cable is not damaged. Check if the connection between the plug and power supply, as well as the cord with the device is correct. Check whether the electrical outlet or car lighter socket are not damaged - check with another source.
If the problem persists, contact the manufacturer's service department.	

11. STORAGE AND TRANSPORT

STORAGE OF EQUIPMENT

If you are not going to use the device for a long time, after disconnecting it from the power supply, cooling and cleaning its components, store it in a dry room away from heat and moisture sources. It is best to place the device in its original packaging, out of reach of children and pets.

TRANSPORT OF THE DEVICE

It is recommended to transport the device in its original packaging. Do not expose the unit to excessive vibration or shock during transportation. Do not move the object when in use. Do it only once it has been switched off and cooled down. Do not expose the device to high temperatures or moisture.

12. WARRANTY AND SERVICE

The product is under warranty for 3 years from the date of sale of the product, except for accessories (AC power supply, car power cord) which are under warranty for 6 months. If you need to make use of a warranty claim, please contact the manufacturer's helpline or point of sale. The product returned for warranty claim should be complete and best in its original packaging.

13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS

mPTech Sp. z o.o. hereby declares that the product complies with all relevant provisions of the so called the 'New Approach' Directive of the European Union.

14. MANUFACTURER

mPTech Sp. z o.o.

Nowogrodzka 31 street
00-511 Warsaw
Poland

Helpline: (+48 71) 71 77 400
E-mail: pomoc@mpotech.eu
Manufactured in China
Lot number: 202110