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Lot number: 612744

**Equipment delivery location:**

mPTech Sp. z o.o.
Krakowska street 119
50-428 Wrocław
Poland

Designed in Poland
Manufactured in China

Manufacturer:

mPTech Sp. z o.o.
Nowogrodzka 31 street
00-511 Warsaw
Poland

Helpline: (+48 71) 71 77 400
E-mail: pomoc@mpotech.eu
Website: www.mpotech.eu

Table of Contents

1. Safety of use.....	8
2. Technical Specification.....	8
2.1 Technical data.....	8
2.2 Package content	9
3. Device construction	9
4. The use of the equipment	9
4.1 Charging the battery.....	9
4.2 Turning the device on, off and operating it	9
4.3 Use	10
5. Troubleshooting	11
6. Cleaning of the device.....	11
7. Correct disposal of used equipment and battery	12
8. The use of the manual	12
9. Warranty and service	12
10. Declaration of conformity with European Union Directives	12

1. Safety of use

Read these guidelines carefully. Not following them may be dangerous or illegal.

–INTENDED USE–

The device is intended for facial skin cleansing using water and ultrasounds and only and for private indoor domestic use only and must not be used for commercial purposes. Any other use is improper and dangerous. The manufacturer is not liable for injuries resulting from the use of the device in a manner inconsistent with its intended use. The device is not a medical equipment.

–CHILDREN SAFETY–

Keep the device, its accessories and packaging out of reach of children. Cables and plastic packaging may wrap around the child's neck or be swallowed, resulting in suffocation or choking.

Children must not play with this device (it is not a toy), or clean it without an adult supervision.

–DO NOT RISK–

Do not use the device near flammable fabrics. Do not use the product with a damaged body or accessories. The product is not intended for use on animals. Do not use the device for purposes other than cleaning the skin.

–QUALIFIED SERVICE–

This product may only be repaired by qualified service personnel of the manufacturer or by an authorised service centre. Repairing the device by an unqualified or unauthorized service technician may result in damage to the device and void the warranty. The manufacturer is not responsible for problems caused by unauthorized modifications to the equipment made by the user.

–CHILDREN AND PHYSICALLY OR MENTALLY DISABLED PEOPLE–

This appliance may be used by children from 8 years of age and older and persons with limited physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised or have been instructed in the safe use of the appliance and understand the hazards involved.

–WATER AND OTHER LIQUIDS–

Protect the device from water and other liquids. Avoid using the device in environments with high humidity like kitchens or bathrooms. Water particles may cause moisture to appear in the device, which adversely affects the electronics and may cause damage to the equipment.

–SURROUNDINGS–

The device should not be used in dusty conditions. While charging the battery make sure it cannot be damaged or endanger people or animals. Do not expose the device to sunlight and do not use it or leave it in close proximity to heat-generating devices (heaters, radiators etc.). Plastic parts may deform, which may affect the operation of the device or cause it to break down completely.

–BATTERY–

Avoid exposing the battery to very high/low temperatures (0°C – 45 °C / 32°F – 113°F). Extreme temperatures may influence capacity and battery life. Avoid contact with liquids and metal objects as it may cause partial or complete battery damage. The battery should be used according to its purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Worn-out or impaired battery should be placed in a dedicated container. Excessive charging the battery can cause damage. Therefore do not charge longer than necessary. The battery in the device is not replaceable.

–CONTRAINDICATIONS–

Do not use the massager without a doctor's approval in the following cases: pregnancy, breastfeeding, cancer, presence of metal bone implants in the areas where the device is used, fever, skin diseases, epilepsy, hypertension, heart failure, purulent dermatitis, anti-acne therapy with retinoids.

If you experience pain or discomfort while using the device, stop immediately.

2. Technical Specification

2.1 Technical data

Device for cavitation peeling	
Model	BE-100224
Power	1.5 W
Power supply	5 V – 500 mA
Battery	Li-ion 3,7 V – 250 mAh, 0,93 Wh
Additional	USB type C cable

2.2 Package content

- Device for cavitation peeling
- USB type C cable
- User Manual


! If any of the above items is missing or damaged, please contact your dealer.

! We advise to keep the container in case of any complaints.

You print a version of the user manual with the large font size or download it from the web address below:

<https://mptech.eu/qr/huslog>

3. Device construction

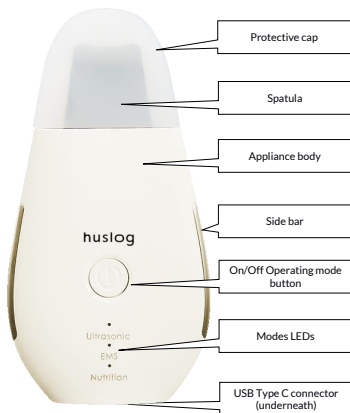4. The use of the equipment

4.1 Charging the battery

! The battery should be fully charged before the first use of product. Several charging and discharging cycles may be necessary to achieve full battery capacity.

! Avoid excessive battery discharge. This can damage the battery or shorten its service life.

! A low battery is indicated by a short double beep and the device is switched off.

! The device can heat up during charging. It is a regular occurrence and should not affect the stability or performance of the device. Do not leave the device to be charged unattended as this may result in fire, damage to the device or its surroundings etc.

To charge your device, connect the USB type C plug of the supplied cable to the USB_C on your device. Connect the other end of the cable to a powered USB adapter or, for example, a socket on your computer, a powerbank, etc., whose output guarantees a power supply of 5 V – 500 mA or higher. When charging starts, the mode LEDs light up and turn off sequentially from the bottom.

When charging is complete all mode LEDs illuminate continuously. Once charged disconnect USB from the device and from the power source.

It is not possible to use the device while the battery is charging.

4.2 Turning the device on, off and operating it

Press and hold the ON/OFF button for approx. 2 seconds to switch the device on or off. Switching on the device will be signalled by a short single beep and all mode LEDs will light up. Further short presses on the ON/OFF button will switch between the three operating modes of the unit:

- **Ultrasonic** - cavitation peeling
- **EMS** - electronic muscle stimulation
- **Nutrition** - the introduction of nutrients deep into the skin

To switch the device off, press and hold the ON/OFF button for approx. 2 seconds. In addition, the product switches off automatically approximately 5 minutes after switching on. Switching off will be signalled by a short single beep.

4.3 Use

 *Avoid using the device for too long and intensively in one place. It can cause skin irritation.*

 *Do not use the device on an irritated or damaged skin.*

 *The device can be used with face, neck and neckline/cleavage.*


 *When using the appliance, always hold it so that it touches the side rails. Do not hold the unit by the body alone.*

Cavitation - a physical phenomenon involving a phase change from a liquid phase to a gas phase. In the case of this product, the cavitation phenomenon occurs when ultrasonic waves are applied to damp skin. This produces microscopic gas bubbles which, when they burst, cause the top layers of the epidermis to peel away along with impurities and excess sebum, revealing a younger layer of skin. Cavitation peeling, when following contraindications (chapter 1. Safety of use), is completely painless and safe.

Clean the skin on your face, neck and neckline before using the device. Remove jewellery. To prevent your hair from falling into your face, secure it with i.e. headband. You can also make a steam for a few minutes or apply a compress with a warm towel - this slowly opens the pores of the skin.

CAVITATION PEELING (ULTRASONIC)


Moisten the skin with warm water or facial cleanser where the spatula of the device is applied. Move the spatula in the direction opposite to the pore alignment, as shown in the illustration below.

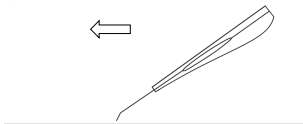While cleansing the nose - both sides and back - move the spatula along the skin upwards. From the jaw to the cheekbone and on the cheeks upwards. On the chin, above the mouth, on the forehead, on the neck and the neckline also upwards. When you have finished exfoliating, wash your skin with clean water or toner. Afterwards, a moisturising cream can be applied to the skin.


ELECTRONIC MUSCLE STIMULATION (EMS)

Moisten your skin with saline or tonic. In this mode of operation, the emission of quick and gentle electrical impulses helps to give the skin elasticity. It also strengthens facial muscles and restores their elasticity. When using the device in this mode, move the spatula according to the direction of the pores, as shown in the illustration below.

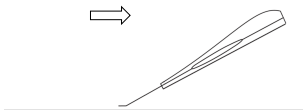NUTRIENTS DEEP INTO THE SKIN (NUTRITION)

Before using the device, apply a vitamin ampoule, serum or cream suitable for this type of treatment. Place the spatula against the skin and move it according to the orientation of the pores, as shown in the illustration below.

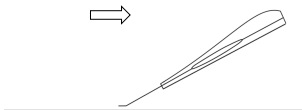 **The skin in the area where the spatula is currently being moved must always be moistened according to the instructions for the mode of operation. Never use the device on dry skin.**

When using the device, do not press the spatula too hard against the skin. Try to keep it tilted at an angle of about 45 degrees to the skin surface. Move the spatula along the skin at a slow pace and do not keep it in one place for too long.

The frequency of application of the device is determined by the individual needs of the skin. We recommend prior consultation with a specialist in this matter. However, the device should not be used several times in one day. Always disinfect the spatula before and after use.

5. Troubleshooting

If you have problems with the unit, please refer to the following tips.

The device does not switch on	Check that the battery is charged. Check that the cables, plugs and sockets supplying power to the unit are not damaged.
After switching on, the spatula does not vibrate	Check that you are not pressing the spatula against the skin too hard. Perhaps there has been a damage to the device. Contact the manufacturer's service department.
If the problem persists, contact the manufacturer's service department.	

6. Cleaning of the device

To extend the service life of the device:

- **WARNING!** The device must not be submerged in water or washed under running water!
- Clean the device after each use.
- Disinfect with any disinfectant or with hydrogen peroxide.
- Before cleaning, make sure that the device is switched off and is not being charged.
- For cleaning use only a dry cloth. Never use any agent with a high concentration of acids, alkalis or abrasives.
- Store the device and its accessories in a dry place away from the reach of children, sources of heat and moisture.
- Avoid very high temperatures, for it can shorten the life of electric and electronic components, melt the plastic parts and destroy the battery.

The packaging of the device can be recycled in accordance with local regulations.

Symbol for corrugated cardboard (packaging material).

The packaging elements should be disposed of in appropriate containers.

Marking of the container into which the stamping should be thrown - plastic/metal.

The marking of the container into which the box - paper should be discarded.

7. Correct disposal of used equipment and battery

The device is marked with a crossed-out garbage bin, in accordance with the European Directive 2012/19/EU on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.

According to the EU directive 2006/66/EC with changes contained in EU directive 2013/56/UE on the disposal of batteries, this product is marked with crossed-out garbage bin symbol. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations. It is not allowed to dispose of accumulators with unsorted municipal waste. Users of battery and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.

8. The use of the manual

Without prior written consent of the manufacturer, no part of this manual may be photocopied, reproduced, translated or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording or storing in any systems of storing and sharing information.

All brand names and product names are trademarks or registered trademarks of their respective companies.

9. Warranty and service

The product is under warranty for 2 years from the date of sale of the product, except for accessories (USB cable) which are under warranty for 6 months. If you need to make use of a warranty claim, please contact the manufacturer's helpline or point of sale. The product returned for warranty claim should be complete and best in its original packaging.

10. Declaration of conformity with European Union Directives

mPTech sp. z o.o. hereby declares that the product complies with all relevant provisions of the so-called 'New Approach' of Directive of the European Union.